Product Manual

ViraSafe™ Lentiviral Bicistronic Expression System (Puro), Pantropic

Catalog Number

VPK-215-PAN 1 kit

FOR RESEARCH USE ONLY Not for use in diagnostic procedures



Introduction

Lentivirus vector based on the human immunodeficiency virus-1 (HIV-1) has become a promising vector for gene transfer studies. The advantageous feature of lentivirus vector is the ability of gene transfer and integration into dividing and non-dividing cells. The pseudotyped envelope with vesicular stomatitis virus envelope G (VSV-G) protein broadens the target cell range. Lentiviral vectors have been shown to deliver genes to neurons, lymphocytes and macrophages, cell types that previous retrovirus vectors could not be used. Lentiviral vectors have also proven to be effective in transducing brain, liver, muscle, and retina *in vivo* without toxicity or immune responses. Recently, the lentivirus system is widely used to integrate siRNA efficiently in a wide variety of cell lines and primary cells both *in vitro* and *in vivo*.

Lentivirus particles are produced from 293T cells through transient transfection of plasmids that encode for the components of the virion (Figure 1). Due to safety concerns regarding the infectious nature of HIV-1, recent lentiviral packaging systems have separated the viral components into 3 or 4 plasmids. However, these systems still present a small chance of generating replication-competent lentivirus upon recombination. In addition, most commercial lentiviral packaging systems provide plasmids containing the viral structure proteins in a premixed formulation, making it nearly impossible to optimize the ratio of the various plasmids for your particular experiment and host cell. Also, most commercial lentivirus transfer vectors contain promoters, antibiotic selection markers and/or reporter genes which may not be optimal or even suitable for your particular expression studies.

Cell Biolabs' ViraSafeTM Lentiviral Expression System provides a much safer method to package lentivirus, while still providing high viral titers. The sequence homology with native HIV-1 has been reduced by 80-90% even compared with other commercial third-generation packaging systems. In addition, each plasmid is provided separately rather than in a packaging mixture. This allows you the flexibility to amplify individual plasmids and optimize the ratio of plasmids for your experiment.

pSMPUW-IRES-Puro Lentiviral Expression Vector contains EF-1 α promoter ahead of the multiple cloning sites, followed by an IRES and puromycin resistant gene (Figure 2).

Key Features of ViraSafe[™] Lentiviral Expression System:

- 1. Transfer Plasmid: Reduce extent of HIV sequences to increase capability up to 10 kb and reduce likelihood of recombination between vector components. Add elements to increase titer and further improve safety.
- 2. Packaging Plasmid: Improve the packaging plasmid to increase performance and reduce the likelihood of recombination between vector components.
 - a. Minimize HIV sequences no accessory proteins, Tat or Rev, or LTRs
 - b. Prevent overlap with vector SM by codon wobbling Gag sequences
 - c. Boost particle production by incorporating adenovirus VA_I element
- 3. Flexible: All vectors including packaging vectors are provided separately to allow end-user to optimize the vector ratio for maximal lentivirus production.



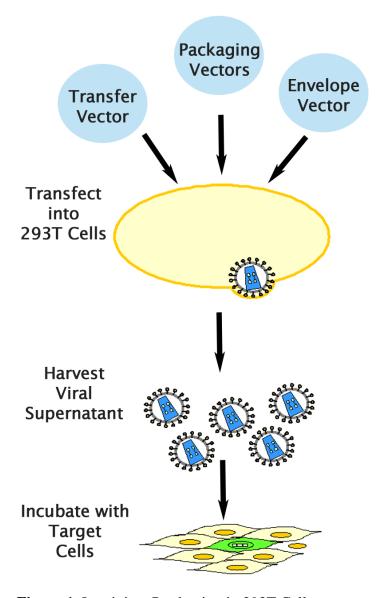


Figure 1. Lentivirus Production in 293T Cells

Related Products

- 1. LTV-100: 293LTV Cell Line
- 2. LTV-200: ViraDuctinTM Lentivirus Transduction Kit
- 3. LTV-300: GFP Lentivirus Control
- 4. VPK-104: ViraBind™ Lentivirus Purification Kit
- 5. VPK-107: QuickTiterTM Lentivirus Titer Kit (Lentivirus-Associated HIV p24)
- 6. VPK-108-H: QuickTiter™ Lentivirus Quantitation Kit (HIV p24 ELISA)
- 7. VPK-205: ViraSafeTM Lentivirus Packaging System, Ecotropic
- 8. VPK-211: pSMPUW Universal Lentiviral Expression Vector (Promoterless)



Unique Elements of the ViraSafeTM Lentivirus Expression System

Third Generation	ViraSafe™	
Lentivirus Expression System	Lentivirus Expression System	
TRANSFER VECTOR		
CMV R U5 SD Agag P MCS RRE R U5	CMV R U5 P MCS WPRE R U5	
PACKAGING VECTOR #1		
P gag RRE Poly(A)	P gag RRE VA Poly(A)	
PACKAGING VECTOR #2		
P Rev	P Rev	
ENVELOPE VECTOR		
P VSVG	P VSVG	

Vector	Element	Name	Benefits compared to other 3 rd Generation Systems
Name			
			ELEMENTS ADDED
Transfer Vector	сРРТ	Central Polypurine Tract	 Increased gene expression levels
		Hybrid 3' LTR Poly(A)	 Increased safety: prevents read-through transcription Increased viral titer: vector transcript more stable in packaging cells
	WPRE	WPRE	Increased viral titer
Packaging Vector #1		Codon Wobble	Increased safety: reduces sequence homology
	VA	Adenovirus VA	Increased viral titer
			ELEMENTS REMOVED
Transfer Vector	Δgag	Gag sequence	Increased safety: reduces sequence homology
	RRE	Rev-Responsive Element	Increased safety: reduces sequence homology

Kit Components

- 1. <u>pSMPUW-IRES-Puro Lentiviral Expression Vector</u> (Part No. VPK-215): One 40 μ L vial at 0.25 mg/mL. The plasmid is kanamycin resistant.
 - Note: Please see Figure 2 for important instructions on bacterial culture of this plasmid.
- 2. pRSV-Rev Packaging Vector (Part No. 320022): One 40 µL vial at 0.25 mg/mL.
- 3. pCMV-VSV-G Envelope Vector (Part No. RV-110): One 40 µL vial at 0.25 mg/mL.
- 4. pCgpV Packaging Vector (Part No. 320024): One 40 μL vial at 0.25 mg/mL.
- 5. <u>pSMPUW-LacZ Control Vector</u> (Part No. 320025): One 40 μL vial at 0.25 mg/mL containing a nuclear localized LacZ driven by MND retroviral LTR promoter. The plasmid is kanamycin resistant.

Note: Please see Figure 2 for important instructions on bacterial culture of this plasmid.

Materials Not Supplied

- 1. 293T cells: we recommend 293LTV Cell Line (Cat.# LTV-100) for high titer production of lentivirus.
- 2. Cell Culture Medium
- 3. Transfection Reagents

Storage

Upon receipt, store all other kit components at -20°C until their expiration dates.

Safety Considerations

Remember that you will be working with samples containing infectious virus. Follow the recommended NIH guidelines for all materials containing BSL-2 organisms. The ViraSafeTM Lentiviral Expression System is designed to minimize the chance of generating replication-competent lentivirus, but precautions should still be taken to avoid direct contact with viral supernatants.

pSMPUW-IRES-Puro Vector

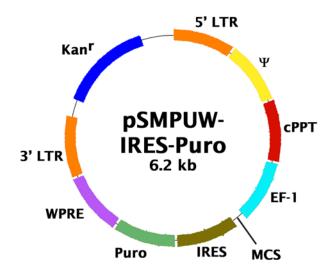


Figure 2: pSMPUW-IRES-Puro Lentiviral Expression Vector (6237 bp, **Kanamycin**-resistant). EcoRI and XhoI Digestion: 1687 bp + 4550 bp.

Note: Bacterial culture of pSMPUW vectors should be done in medium containing 10 µg/mL Kanamycin. For maximal plasmid yield and quality, we recommend Stbl3 endoA1+ competent cells (Invitrogen) and treatment with alkaline proteinase (Promega #A1441 or Sigma #P8038) for 4-5 min using 10 units of proteinase per mL of bacterial lysate before adding neutralization solution.

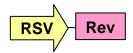


Figure 3: pRSV-Rev Packaging Vector (4180 bp, **Ampicillin**-resistant). EcoRI Digestion: 300 bp + 3880 bp

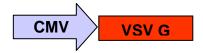


Figure 4: pCMV-VSV-G Envelop Vector (6051 bp, **Ampicillin**-resistant). EcoRI Digestion: 787 bp + 1668 bp + 3596 bp.

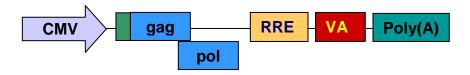


Figure 5: pCgpV Packaging Vector (9118 bp, **Ampicillin**-resistant). Pst I Digestion: 927 bp + 1424 bp + 6767 bp.



Lentivirus Production

- 1. One day before transfection, plate sufficient 293T cells or 293LTV cells (Cat. # LTV-100) to achieve 70-80% confluence on the day of transfection.
- 2. Transfect cells by Calcium Phosphate or other transfection reagents.

Note: We suggest transfecting cells with FuGENE® Transfection Reagent (Roche Applied Science) or LipofectamineTM Plus (Invitrogen). We recommend the ratio of vectors at 3:1:1:1 (pSMPUW: pCMV-VSV-G:pRSV-REV:pCgpV).

- 3. Harvest lentiviral supernatant 36-72 hours after transfection. Supernatant can be harvested 2 or 3 times, every 12 hours. Keep it at 4°C over the collecting period.
- 4. Pool the collected supernatants, centrifuge 5 minutes at 1500 rpm to remove cell debris and filtrate on 0.22 μm.
- 5. Supernatants can be used directly or purified/concentrated if needed. For long term storage, store supernatant at -80°C in aliquots.

Post-Packaging Considerations

Packaging your lentivirus is only the first step to ensuring successful expression of your gene. The following steps should be considered prior to infection of your host cell:

- 1. **Concentration and purification of your lentivirus**: Because of the latent nature of lentivirus, it is imperative that your virus be highly concentrated before infecting your host cell. Also, impurities from your viral supernatant can decrease the efficiency of infection. We recommend using Cell Biolabs' ViraBindTM Lentivirus Concentration and Purification Kit (Catalog # VPK-090).
- 2. **Measure the titer of your lentivirus**: This is an important step to ensure consistent viral transduction into your host cell. However, QPCR or stable clone counting can take as much as 1-2 weeks to perform. Traditional p24 ELISA kits can greatly overestimate your lentiviral titer. Our advanced p24 ELISA, QuickTiterTM Lentivirus Titer Kit (Catalog # VPK-107), uses exclusive technology that eliminates free p24 from your supernatant, giving you much more accurate lentiviral titers. Results are obtained in 6-18 hours.
- 3. Use transduction reagents to increase infection efficiency: Many cells are difficult to infect with lentivirus, and without supplemental reagents transduction efficiencies can be low. Reagents such as Polybrene® can help, but are often insufficient. Cell Biolabs' proprietary reagents in our ViraDuctinTM Lentivirus Transduction Kit (Catalog # LTV-200) form a super-complex with your virus to increase transduction efficiencies by promoting virus and cell interaction.



Example of Results

The following figure demonstrates typical results seen with Cell Biolabs ViraSafeTM Lentiviral Expression System. One should use the data below for reference only.

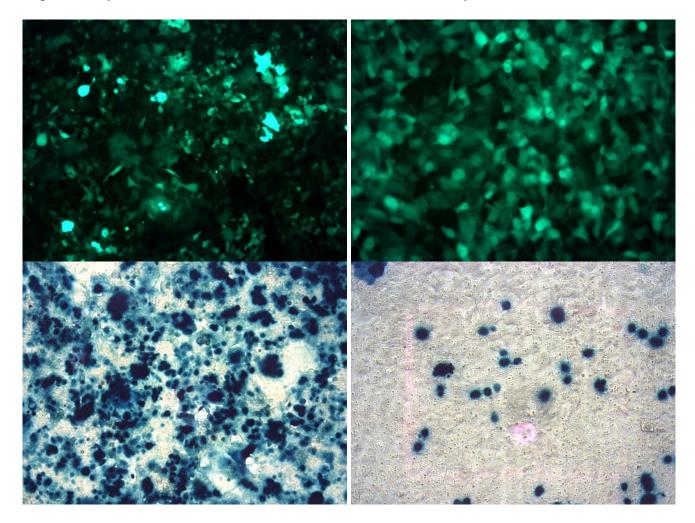


Figure 6: GFP and nLacZ Lentivirus Production and Transduction: Lentiviral supernatant is produced by cotransfecting 293T cells (Cat.# LTV-100) with pLenti-GFP (Cat.# LTV-400) or pSMPUW-MNDnLacZ (Cat.# LTV-402) and ViraSafeTM Lentiviral Packaging System (Cat. # VPK-206). 293AD cells (Cat.# AD-100) are seeded at 100,000 cells/well in a 6-well plate overnight. Cells were infected with GFP or nLacZ lentivirus in the presence of 8 μg/mL Polybrene for 72 hrs. **Left:** 293LTV Transfection; **Right:** 293AD Transduction.

Appendix

pSMPUW-IRES-Puro Plasmid Sequence

Pink: 5' CMV/LTR, ψ, cPPT

Blue: EF-1
Purple: MCS
Red: IRES
Green: Puro
Brown: WPRE
Orange: 3' LTR

Blue: Kanamycin Resistance gene

CAAGTGTATCATATGCCAAGTACGCCCCCTATTGACGTCAATGACGGTAAATGGCCCGCCTGGCATTATGCCCAGTACATGACCTTATGGGACTTTCCTACTTGGGGAAAATAGCGGCCGCCACAATTTTAAAAGAAAAGGGGGGGATTGGGGGGTACAGTGCAGGGGAAAGAATAGTAGACATAATAGCAACAGACATACAAACTA ACAAAAGGACTCGCCCTGCCTTGGGGAATCCCAGGGACCGTCGTTAAACTCCCACTAACGTAGAACCCAGAGATCGCTGCGTTCCCGCCCCTCACCCGCCCGCTCTCGTCATCACTGAGGTGGAGAAGAGCATGCGTGAGGCTCCGGTGCCCGTCAGTGGCCAGAGCGCACATCGCCCACAGTCCCCGAGAAGTTGGGGGGGAGG ${\tt GCCGAAGCCGCTTGGAATAAGGCCGGTGTGCGTTTGTCTATATGTTATTTTCCACCATATTGCCGTCTTTTGGCAATGTGAGGGCCCGGAAACCTGGCCCTGT}$ ACCCCATTGTATGGGGATCTGATCTGGGGCCTCGGTGCACATGCTTTACATGTGTTTAGTCGAGGTTAAAAAAACGTCTAGGCCCCCCGAACCACGGGGACGTG GTTTTCCTTTGAAAAACACGATGATAATATGGCCACAACCATGGTTACCGAGTACAAGCCCACGGTGCGCCTCGCCACCCGCGACGACGTCCCCAGGGCCGTA TGTGGATACGCTGCTTTAATGCCTTTGTATCATGCTATTGCTTCCCGTATGGCTTTCATTTTCTCCTCCTTGTATAAATCCTGGTTGCTGTCTCTTTATGAGGAGTTGTGGCCCGTTGTCAGGCAACGTGGCGTGGTGTGCACTGTGTTTGCTGACGCAACCCCCACTGGTTGGGGCATTGCCACCACCTGTCAGCTCCTTTCCGG $\tt CAGCGGACCTTCCTTCCCGCGGCCTGCTGCCGGCTCTTCCGCGTCTTCCGCGTCTTCGCCTCAGACGACTGGATCTCCCTTTGGGCCGCCTCCCC$ GCTTAGTACTGGTACCTTTAAGACCAATGACTTACAAGGCAGCTGTAGATCTTAGCCACTTTTTAAAAGAAAAGGGGGGACTGGAAGGGCTAATTCACTCCCA AGTAGTGAGGAGGCTTTTTTTGGAGGCCTAGGCTAGAGATCATAATCAGCCATACCACATTTGTAGAGGTTTTACTTGCTTTAAAAAACCTCCCACACCTCCCC $\tt CTGAACCTGAAACATAAAATGAATGCAATTGTTGTTGTTTAACTTGTTTATTGCAGCTTATAATGGTTACAAATAAAGCAATAGCATCACAAATTTCACAAATA$ $\tt CTGCCGCTTACCGGATACCTGTCCGCCTTTCTCCCTTCGGGAAGCGTGGCGCTTTCTCATAGCTCACGCTGTAGGTATCTCAGTTCGGTGTAGGTCGTTCGCT$



References

- 1. Chen, M. et al. (2002). *Nature Genetics* **32(4)**: 670-675.
- 2. Naldini, L., U. Blomer, P. Gallay, D. Ory, R. Mulligan, F. H. Gage, I. M. Verma, and D. Trono (1996) *Science* **272**:263-267.
- 3. Verma, I. M., and N. Somia (1997) *Nature* **389**:239-242
- 4. Kahl C. A., Marsh J., Fyffe J., Sanders D. A., and K. Cornetta (2004) *J Virol.* **78**:1421-30.
- 5. White S. M., Renda M., Nam N. Y., Klimatcheva E., Zhu Y., Fisk J., Halterman M., Rimel B. J., Federoff H., Pandya S., Rosenblatt J. D., and V. Planelles (1999) *J Virol.* **73**:2832-40.
- 6. Kafri T., van Praag H., Ouyang L., Gage F. H., and I. M. Verma (1999) *J Virol.* **73**:576-84.

Notice to Purchaser

This product is sold for research and development purposes only and is not to be incorporated into products for resale without written permission from Cell Biolabs. The patented technology is covered by a license from CHLA and University of Southern California. By the use of this product you accept the terms and conditions of all applicable Limited Use Label Licenses. You may contact our Business Development department at busdev@cellbiolabs.com for information on sublicensing this technology.

Warranty

These products are warranted to perform as described in their labeling and in Cell Biolabs literature when used in accordance with their instructions. THERE ARE NO WARRANTIES THAT EXTEND BEYOND THIS EXPRESSED WARRANTY AND CELL BIOLABS DISCLAIMS ANY IMPLIED WARRANTY OF MERCHANTABILITY OR WARRANTY OF FITNESS FOR PARTICULAR PURPOSE. CELL BIOLABS' sole obligation and purchaser's exclusive remedy for breach of this warranty shall be, at the option of CELL BIOLABS, to repair or replace the products. In no event shall CELL BIOLABS be liable for any proximate, incidental or consequential damages in connection with the products.



Contact Information

Cell Biolabs, Inc. 7758 Arjons Drive San Diego, CA 92126

Worldwide: +1 858-271-6500 USA Toll-Free: 1-888-CBL-0505 E-mail: tech@cellbiolabs.com

www.cellbiolabs.com

©2009-2015: Cell Biolabs, Inc. - All rights reserved. No part of these works may be reproduced in any form without permissions in writing.

